



iQuant™ 1X dsDNA BR Assay Kit (2 – 1000 ng)

Catalog Number: N026, N027

Table 1. Kit Components and Storage

Material	Amount	Concentration	Storage	Stability
iQuant™ 1X dsDNA BR Assay Kit (Cat. No. N026, 200 assays)				The product is stable for at least 6 months when stored as directed.
iQuant™ 1X dsDNA BR Reagent (Component A)	40 mL	1X	2-8 °C Protect from light	
dsDNA Standard #1 (Component B)	0.4 mL	0 ng/μL in TE buffer		
dsDNA Standard #2 (Component C)	0.4 mL	100 ng/μL in TE buffer		
iQuant™ 1X dsDNA BR Assay Kit (Cat. No. N027, 500 assays)				
iQuant™ 1X dsDNA BR Reagent (Component A)	100 mL	1X	2-8 °C Protect from light	
dsDNA Standard #1 (Component B)	1 mL	0 ng/μL in TE buffer		
dsDNA Standard #2 (Component C)	1 mL	100 ng/μL in TE buffer		

Approximate fluorescence excitation/emission maxima, in nm: 515/535, bound to DNA.

Product Description

The iQuant™ 1X dsDNA BR Assay Kit provides a simple, sensitive, and accurate quantitation for dsDNA. The kit contains ready-to-use assay reagent, and pre-diluted dsDNA standards. The assay kit is highly selective for dsDNA, and highly reliable in detecting dsDNA ranging from 2 to 1000 ng, and offers advantages in stability, linear dynamic range, and sensitivity over other traditional of DNA quantitation. The assay is performed at room temperature. Simply add your sample (any volume between 1 μl and 20 μl is acceptable) to the assay reagent, and read the fluorescence using fluorescence plate reader or Fluorometer such as Qubit® Fluorometer. The kit is well tolerated to common contaminants such as proteins, salts, solvents and detergents.

Handling and Disposal

There is no safety data available for iQuant™ 1X dsDNA BR reagent. Treat the iQuant™ 1X dsDNA BR reagent with the safety precautions as other potentially harmful reagents and to dispose of the reagent in accordance with local regulations. Centrifuge the dsDNA standards before opening vials to minimize loss on the cap. Use properly calibrated pipettes for best accuracy.

General Protocol

1. Measure dsDNA samples using a Fluorescence Microplate Reader

(Note: For simplicity, the following protocol is written using 10 μL of dsDNA sample volume. In practice, the volume of dsDNA sample could be ranging from 1 μL to 20 μL depending on the concentration of dsDNA sample, then adjust the volume of iQuant™ 1X dsDNA BR reagent to 200 μL.)

1.1 Warm up the iQuant™ 1X dsDNA BR Assay Kit to room temperature.

- 1.2 Add 190 μL of the iQuant™ 1X dsDNA BR reagent (Component A) to each well of a black 96-well microplate. Black plates such as Greiner or Corning black 96-well plates are recommended to minimize fluorescence bleed-through from other well.
- 1.3 Prepare a series of dsDNA standard dilutes from dsDNA Standard #2 (Component C) or your known dsDNA sample.
- 1.4 Add 10 μL of each dsDNA standard dilutes and the unknown dsDNA samples in duplicate or triplicates into separated wells and mix well by pipetting up and down.
- 1.5 Incubate the microplate at room temperature for 2 minutes in the dark.
- 1.6 Measure the fluorescence using a microplate reader with 485 nm excitation and 530 nm emission, with the appropriate cut-off.
- 1.7 Generate a linear standard curve by plotting fluorescence versus DNA concentration of the DNA standards. Use the standard curve and the fluorescence of the unknown DNA samples to determine the unknown DNA concentration.

2. Measure dsDNA samples using the Qubit® Fluorometer from ThermoFisher

(Note: For simplicity, the following protocol is written using 10 μL of dsDNA sample volume. In practice, the volume of dsDNA sample could be ranging from 1 μL to 20 μL depending on the concentration of dsDNA sample, then adjust the volume of iQuant™ 1X dsDNA BR reagent to 200 μL .)

- 2.1. Warm up the iQuant™ 1X dsDNA BR Assay Kit to room temperature.
- 2.2. Add 190 μL of the iQuant™ 1X dsDNA BR reagent (Component A) to each assay tube. (**Note:** Use only thin-wall, clear 0.5 mL PCR tubes. iQuant™ assay tubes (Cat No. N022)).
- 2.3. Add 10 μL of dsDNA standard #1 (Component B), dsDNA standard #2 (Component C), and the unknown dsDNA samples to the appropriate tubes and mix by vortexing 2-3 seconds, and label the lids of each DNA standard tube and unknown sample tubes correctly.
- 2.4. Incubate all tubes at room temperature for 2 minutes in the dark.
- 2.5. Measure the fluorescence on the Qubit® fluorometer using the **1X dsDNA Broad Range** program, according to the manufacture's recommendation.

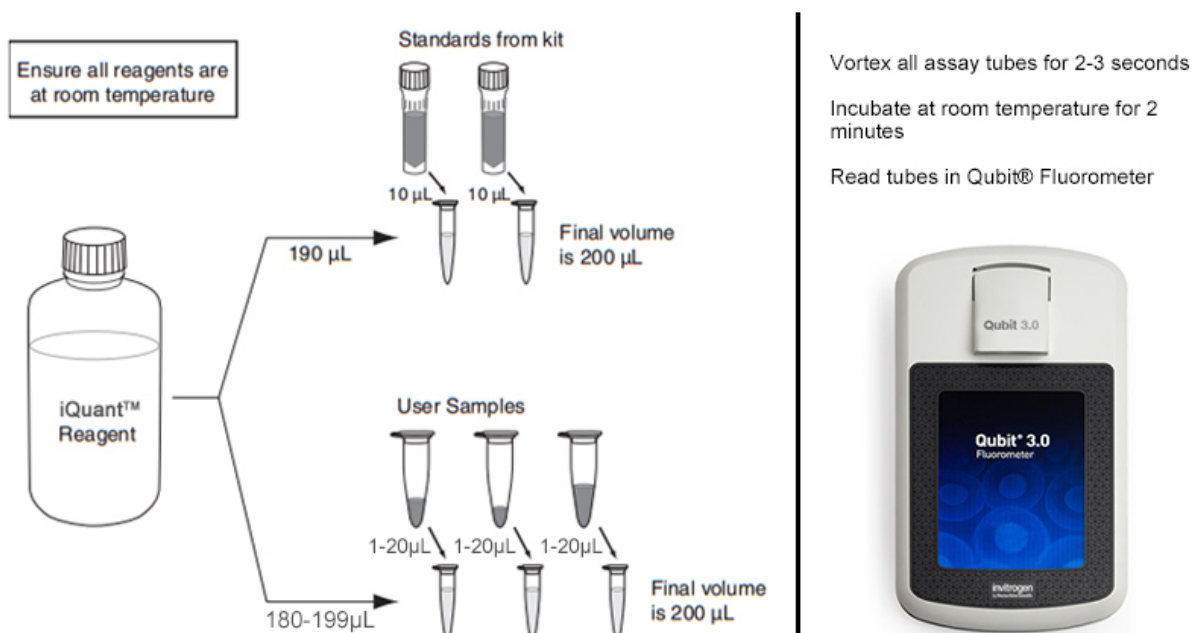


Figure 1. Qubit® assay workflow

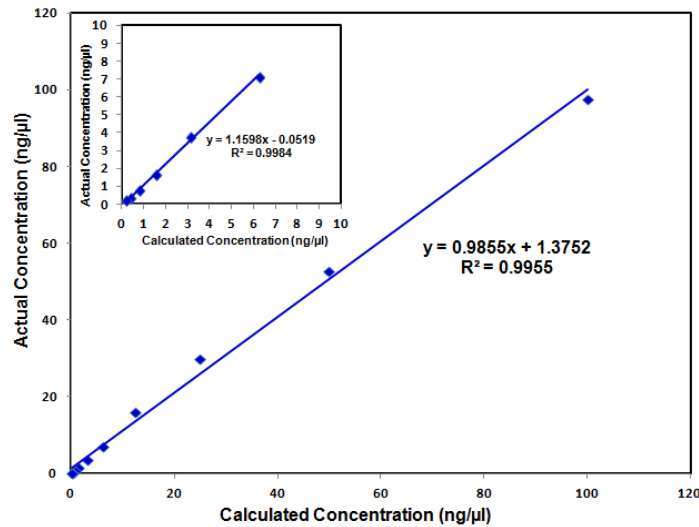


Figure 2. Quantitation of dsDNA using Qubit® Fluorometer.

Considerations for Data Analysis

It is more prefer to use a dsDNA standard similar to the unknown samples (i.e. similar in size, linear vs circular). We found using the iQuant™ 1X dsDNA BR reagent most linear dsDNA yield similar results. If the fluorescence of an unknown sample is higher than dsDNA standard #2 (Component C), further dilute the sample and add 1~10 μL of diluted sample to perform the assay.

Appendix

Table 2. Effect of Contaminants in the iQuant™ 1X dsDNA BR Assay

Contaminant	Final Concentration in Assay	Concentration in 10 μL Sample	Result
Proteins			
Bovine Serum Albumin	10 mg/mL	200 mg/mL	OK
Salts			
Sodium Chloride	20 mM	400 mM	OK
Magnesium Chloride	5 mM	100 mM	OK
Sodium Acetate	20 mM	400 mM	OK
Ammonium Acetate	20 mM	400 mM	OK
Organic Solvents			
Ethanol	0.5%	10%	OK
Chloroform	0.5%	10%	OK
Phenol	0.1%	2%	OK
Detergents			
Sodium Dodecyl Sulfate	0.01%	0.2%	OK
Triton X-100	0.01%	0.2%	OK
Other Compounds			
dNTPs	100 μM	2 mM	OK
RNA	1X	1X	OK
Polyethylene Glycol	1%	20%	OK
Agarose	0.1%	2%	OK

Table 3. iQuant™ assay kits compatible with the Qubit® Fluorometer

Product	Cat. No.	Unit	Target	Notes
iQuant™ 1X dsDNA HS Assay Kit	N020	200T	dsDNA	Detection range: 10 pg/μL – 100 ng/μL; Useful for quantitation of PCR products, viral DNA, and samples for subcloning.
	N021	500T		
iQuant™ dsDNA HS Quantitation Kit	N010	200T	dsDNA	Detection range: 10 pg/μL – 100 ng/μL; Useful for quantitation of PCR products, viral DNA, and samples for subcloning.
	N011	1000T		
iQuant™ 1X dsDNA BR Assay Kit	N026	200T	dsDNA	Detection range: 100 pg/μL – 1000 ng/μL; Useful for quantitation of genomic and miniprep DNA samples.
	N027	500T		
iQuant™ dsDNA BR Quantitation Kit	N012	200T	dsDNA	Detection range: 100 pg/μL – 1000 ng/μL; Useful for quantitation of genomic and miniprep DNA samples.
	N013	1000T		
iQuant™ ssDNA Quantitation Kit	N014	200T	ssDNA	Detection range: 50 pg/μL – 200 ng/μL; Useful for quantitation of oligos, primers, denatured DNA, PCR products.
	N015	1000T		
iQuant™ RNA HS Quantitation Kit	N016	200T	RNA	Detection range: 250 pg/μL – 100 ng/μL; Useful for quantitation of samples of microarray and RT-PCR.
	N017	1000T		
iQuant™ RNA BR Quantitation Kit	N018	200T	RNA	Detection range: 1 ng/μL – 1000 ng/μL; Useful for quantitation of samples of microarray and RT-PCR.
	N019	1000T		

Related Products

Cat. No.	Product Name	Unit Size
N010	iQuant™ dsDNA HS Assay Kit	200 assays
N011	iQuant™ dsDNA HS Assay Kit	1000 assays
N012	iQuant™ dsDNA BR Assay Kit	200 assays
N013	iQuant™ dsDNA BR Assay Kit	1000 assays
N014	iQuant™ ssDNA Assay Kit	200 assays
N015	iQuant™ ssDNA Assay Kit	1000 assays
N016	iQuant™ RNA HS Assay Kit	200 assays
N017	iQuant™ RNA HS Assay Kit	1000 assays
N018	iQuant™ RNA BR Assay Kit	200 assays
N019	iQuant™ RNA BR Assay Kit	1000 assays
N020	iQuant™ 1X dsDNA HS Assay Kit	200 assays
N021	iQuant™ 1X dsDNA HS Assay Kit	500 assays
N026	iQuant™ 1X dsDNA BR Assay Kit	200 assays
N027	iQuant™ 1X dsDNA BR Assay Kit	500 assays
N022	iQuant™ Assay Tubes	500 tubes
N023	iQuant™ microRNA Assay Kit	200 assays
N024	iQuant™ microRNA Assay Kit	1000 assays
N025	iQuant™ RiboGreen RNA Assay Kit	1000 assays